



Orthopaedic Devices, Foot & Ankle Products

Preventive, Postoperative, Off-loading, Wound Care



DARCO



DARCO – Your partner for footcare products

DARCO is an international group operating worldwide in a leading position with orthopaedic cast boots, walkers, postoperative shoes and off-loading boots.

DARCO has become a specialist in diabetic and wound care due to many years of experience. For more than 30 years we represent innovative products.

DARCO products are readily available, offer great value and can be effectively placed on patients and/or modified in minutes. With the rise of diabetes escalating worldwide, the incidence of wounds and ulcerations are increasing as well.

In this respect our R&D team consisting of orthopaedic technicians and product designers focuses on the development of new products representing state-of-art technology in close cooperation with leading medical specialists and health care professionals all over the world.

DARCO is specifically involved in promoting and expanding existing knowledge as well as organizing multiple further training opportunities which include workshops, congresses, symposiums and exhibitions. Over the last years this has become a source of information for medical specialists and a platform for mutual exchange between the USA, Germany, Great Britain and China. With the inauguration of our subsidiary in Dharwad, Karnataka, India became part of this global platform.

This is the way we include practical know-how and modern technology in our product design and how we remain standing firmly on our own two feet.



By courtesy of Dr. Sindhya S. Sukkarni M.D.
Chetana Diabetes Foot Care Center



4 Indications

Postoperative Care

- | | |
|-------------------|---------------|
| 7 MedSurg™ | 9 OrthoWedge® |
| Original MedSurg™ | 10 TwinShoe |
| 8 Softie™ Shoe | 10 TAS™ |
| APB™ | 11 GTS™ |
| 9 HeelWedge™ | |

Trauma Care

- | | |
|--------------------------|--------------------------|
| 13 Body Armor® Cast Shoe | 14 AirTraveler™ |
| SlimLine™ | 15 Body Armor® Walker II |
| 14 FXPro™ Walker | 15 Body Armor® ProTerm |

Off-loading and Wound Care

- | | | |
|----------------------|----------------|-----------------------|
| 17 PegAssist™ Insole | 19 OrthoWedge® | 21 GentleStep™ Insole |
| WCS™ | TwinShoe | WCS™ /-Light Insoles |
| 18 WCS™ Closed-Toe | 20 MedSurg™ | |
| HeelWedge™ | GentleStep™ | |

Plantar Fasciitis Products

- | | |
|---------------------------|----------------------------------|
| 23 Pneumatic Night Splint | 24 DCS™ Plantar Fasciitis Sleeve |
| Body Armor® Night Splint | ProMotion® Plus |

Ankle Supports

- 26 Body Armor® Vario
Web™ Ankle Brace

Pediatric Products



- | | |
|-----------------------------|------------------------|
| 28 Pediatric MedSurg™ | 28 Pediatric SlimLine™ |
| Pediatric Original MedSurg™ | |

- 30 Sizes and Measurements
31 Sizes and Measurements Kids
32 Article List

Indications

		Postoperative Healing	Trauma Care	Fractures	Sprains	Weak Ankles	Plantar Fasciitis	Wounds / Off-loading	Presence of K-Wires	Cast Protection	Diabetic Foot / Neuropathy	All Weather Protection	Arch Support
Postoperative Care													
MedSurg™ (p. 7)		✓	✓						✓	✓	✓		
Original MedSurg™ (p. 7)		✓	✓							✓	✓		
Softie™ Shoe (p. 8)		✓	✓							✓	✓		
APB™ (p. 8)		✓	✓						✓	✓	✓	✓	✓
HeelWedge™ (p. 9)		✓						✓					
OrthoWedge® (p. 9)		✓	✓					✓	✓				
TAS™ (p. 10)		✓											
GTS™ (p. 11)		✓											
Trauma Care													
Body Armor® Cast Shoe (p. 13)				✓							✓	✓	✓
SlimLine™ (p. 13)			✓	✓						✓	✓	✓	
FXPro™ Walker (p. 14)		✓	✓	✓	✓	✓							
AirTraveler™ (p. 14)		✓	✓	✓	✓	✓							
Body Armor® Walker II (p. 15)		✓	✓	✓	✓	✓		✓					
Body Armor® ProTerm (p. 15)		✓	✓										
Off-loading & Wound Care													
PegAssist™ Insole (p. 17)		✓	✓					✓					
WCS™ (p. 17)			✓					✓				✓	

		Postoperative Healing	Trauma Care	Fractures	Sprains	Weak Ankles	Plantar Fasciitis	Wounds / Off-loading	Presence of K-Wires	Cast Protection	Diabetic Foot / Neuropathy	All Weather Protection	Arch Support
WCS™ Closed-Toe (p. 18)			✓						✓			✓	
HeelWedge™ (p. 18)		✓							✓				
OrthoWedge® (p. 19)		✓	✓						✓	✓			
MedSurg™ (p. 20)		✓	✓						✓	✓		✓	
Gentle Step™ (p. 20)												✓	
Gentle Step™ Insole (p. 21)												✓	
WCS™ / -Light Insoles (p. 21)			✓						✓			✓	
Plantar Fasciitis Products													
Pneumatic Night Splint (p. 23)								✓					
Body Armor® Night Splint (p. 23)								✓					
DCS™ Plantar Fasciitis Sleeve (p. 24)								✓					✓
ProMotion® Plus (p. 24)								✓					✓
Ankle Supports													
Body Armor® Vario (p. 26)		✓	✓			✓	✓						✓
Web™ Ankle Brace (p. 26)						✓	✓						
Pediatric Products													
Pediatric MedSurg™ (p. 28)		✓	✓						✓	✓		✓	
Pediatric Original MedSurg™ (p. 28)		✓	✓							✓		✓	
Pediatric SlimLine™ (p. 28)			✓	✓						✓	✓	✓	



Postoperative Care

MedSurg™

Postoperative Shoe




compatible with:

PegAssist™ Insole
for selective off-loading



exchangeable
insole

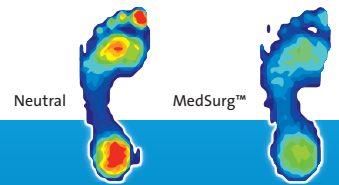


> used postoperatively for stabilization and off-loading

- > 3 Women-Sizes: WS, WM, WL
- > 4 Men-Sizes: MS, MM, ML, MXL
- > right and left fit

The MedSurg™ Shoe with its famous semi-rigid sole is built around our patented MetaShank™ forefoot support, that provides more rigid control and better protection after osteotomies. Our breathable, mesh upper conforms to bony abnormalities of the foot and has a reinforced heel counter and padded heel collar to stabilize the rear foot and reduce slippage. Extra-long straps accommodate even the bulkiest dressings. The MedSurg™ shoe is available in 7 different sizes. Helping you get the best possible fitting post-operative shoe.

Original MedSurg™

Postoperative Shoe

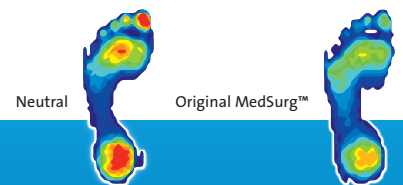

> used postoperatively for stabilization and off-loading

- > 3 Women-Sizes: WS, WM, WL
- > 4 Men-Sizes: MS, MM, ML, MXL
- > right and left fit

The DARCO medical-surgical shoe is the shoe that started it all. Our Original MedSurg™ shoe was invented over 25 years ago and has since become known as the industry standard. Often referred to generically as »the DARCO shoe«, it is anything but generic and its long list of features and excellent quality make it one of the best values in the business.

Softie™ Shoe

Postoperative Shoe



> light and soft postoperative shoe for soft tissue procedures

> 3 Women-Sizes: S, M, L


> 4 Men-Sizes: S, M, L, XL

> right and left fit

The DARCO Softie™ Shoe was designed specifically for use after soft tissue procedures or trauma to the foot. This shoe is a very good solution for those looking for a comfortable, low cost, postoperative shoe. The Semi-rigid metatarsal shank provides superior control of forefoot motion to ensure postoperative healing undisturbed. The Tri-laminated upper is made of light weight breathable tricot-foam and comfortably conforms to bony abnormalities. Hook and Loop Straps are easy to fasten and unfasten with one hand.

APB™

All Purpose Boot

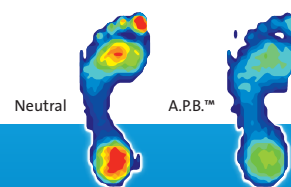

> after surgery and for stabilizing
> it can also serve as a cast boot for particular applications
> impact-absorbing effect and pressure redistribution

> 5 Sizes: XS, S, M, L, XL

> right and left fit

In the A.P.B.™ All Purpose Boot the same sole technology is used as in the proven MedSurg post surgical shoe, so that its use as a plaster shoe or walking sole in rigid bandages is certainly possible. The new closure technique allows for simple and two-sided application. The shoe is completely closed and water-repellent and keeps dressings dry.

compatible with:

PegAssist™ Insole
for selective off-loading



exchangeable
insole

HeelWedge™

Off-loading Shoe



compatible with:




TwinShoe
for height adjustment



PegAssist™ Insole
for selective off-loading

exchangeable
insole

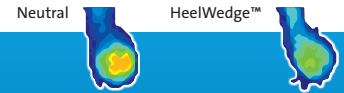

> for the treatment of plantar fasciitis, ulcerations, infections, trauma and following surgery of the soft tissue or bony structure of the heel

- > 5 Sizes: XS, S, M, L, XL
- > right and left fit

The HeelWedge™ is clinically proven to off-load pressure by over 26 % from the heel by shifting weight to the mid and forefoot to promote faster healing after surgery, trauma or when wounds or ulcerations are present on the heel. The square-toe design offers more space for bandage material. The modified and removeable forefoot closure provides better left / right fit. The removable EVA-insole can be modified as needed.

OrthoWedge®

Off-loading Shoe



compatible with:



TwinShoe
for height adjustment



PegAssist™ Insole
for selective off-loading

exchangeable
insole



> use after surgeries to correct hallux valgus, hammer toes /
Tailor's bunions
> treatment of diabetic ulcerations

- > 5 Sizes: XS, S, M, L, XL
- > right and left fit

The OrthoWedge® Off-loading Shoe is specifically designed to protect the forefoot by removing most of the pressure from the metatarsal head and digits during the weight-bearing portion of the gait cycle. For modification the insole can be taken out. The patented construction of the wedge obtains an off-loading mechanism (TPR-sole). The weight is noticeably reduced at the forefoot. Hook and loop fastener traps make opening easy and allow a comfortable fitting.

TwinShoe

Balance Shoe



compatible with:

HeelWedge™, OrthoWedge®,
Body Armor® Walker II

- > for a quick mobilization after foot surgery
- > height adjusting balance shoe matches all orthopaedic devices with a thick sole
- > 5 Sizes: XS, S, M, L, XL
- > right and left fit

The TwinShoe is a matching and height adjusting shoe supplementing our classic product line (listed bottom left) and used specifically for clinical care support. No more time-intensive applying and removing of height adjustments. DARCO TwinShoe warrants a safe and physiologically correct gait.

TAS™

Toe Alignment Splint



- > postoperative after hallux valgus (1), hammer toes (2) and Tailor's bunion (3)
- > 1 Size: Universal
- > right and left fit

The TAS™ Toe Alignment Splint is an effective and inexpensive way of maintaining correct alignment of the toes during the healing process. A thin elastic metatarsal band reduces slippage and soft toe straps hold the toes in the correct position to maintain alignment following surgery to correct hallux valgus, hammer toes or Tailor's bunions. May also be used in the treatment of capsulitis plantar to the m.p. joints caused by mild subluxations of the joint.

Great Toe Splint



> Indicated for relief of symptoms associated with hallux valgus and other great toe joint conditions

- > maintains alignment
- > 24 hour support
- > helps relieve painful bunions
- > can be worn postoperatively
- > can be worn under a sock to accommodate most footwear
- > 2 Women-Sizes: S, L
- > 2 Men-Sizes: S, L

The GTS™ – Great Toe Splint is an adjustable, postoperative splint. Indicated for relief of symptoms associated with hallux valgus and other great toe joint conditions. Manufactured with lightweight neoprene and nylon, the GTS™ can be worn day and night, and its form-fitting design easily accommodates suitable footwear.



Trauma Care

Body Armor® Cast Shoe

Cast Boot



> usable as a walking shoe over a cast bandage

> 5 Sizes: XS, S, M, L, XL

> right and left fit

The new Body Armor® Cast Shoe provides the ultimate cast protection in a new modern design that will appeal to all patients. Soft, durable EVA is combined with a patent protected unique bungee closure system to seat and secure the cast in the shoe, all the while providing the comfort patients need and deserve.

SlimLine™

Cast Boot



> used over a fiberglass- or plaster cast


> shock absorbing effect

> 5 Sizes: XS, S, M, L, XL

> right and left fit

compatible with:

PegAssist™ Insole
for selective off- loading

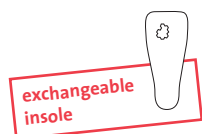

The SlimLine™ is able to enclose voluminous bandages due to its square shape in the forefoot area and a flexible closure technology. For modification the insole can be taken out. The upper is Nylon and the special sole warrants a physiological correct gait.

FXPro™ Walker

Stirrup Walker



- > for tendon sprains and ankle injuries
- > 5 Sizes: XS, S, M, L, XL
- > right and left fit
- > available in high and low versions

The FXPro™ Walker is built around an ultra low-profile sole that provides the smoothest possible transition from heel to toe for normal ambulation without stressing the foot and ankle. Spatula stirrups conform to the calf to reduce chafing and the vented liner ensures that your patients stay cooler and dryer. Circular velcro straps warrant best fit in the leg and are individual adjustable. The rocker sole warrants a physiological correct gait.

AirTraveler™

Pneumatic Stirrup Walker



- > for ruptures and ankle injuries
- > quick and easy with air-pressure
- > 5 Sizes: XS, S, M, L, XL
- > right and left fit
- > available in high and low versions

The AirTraveler™ is an economical, pneumatic walker. In comparison with the FXPro™ Walker, the AirTraveler™ ensures more compression with the circumferential air bladders. Posterior flow restrictor channels forward the air for the greatest immobilizing effect.

Body Armor® Walker II

Cast Replacement



compatible with:

TwinShoe

for height adjustment



> total contact system instead of a cast

> 5 Sizes: XS, S, M, L, XL

> right and left fit

> available in high and low versions

With the hard shell circumferential protection the Body Armor® Walker II can be used like a cast. Leverage fastener similar to ski boot buckles are attached directly at the shell and controls compression and immobilisation. Two insoles can be taken out for further rework and off-loading. An additional (protection) plate at the heel increases the stabilization of the foot enormously.

Body Armor® ProTerm

Foot Stump Orthosis



> interim treatment to support normal stump formation

> rear weight-bearing wedge-shaped outer and innersole to support dorsal flexion

> used postoperatively following Chopart and Lisfranc amputations

> wide opening for easy foot insertion

> practical & secure closure system

> air compression bladder for customizable fitting

> 5 Sizes: XS, S, M, L, XL

> 2 Sizes: M, L

> can be worn on the right and left



Off-loading and Wound Care


PegAssist™ Insole

Off-loading Insole



compatible with:

MedSurg™, A.P.B.™, HeelWedge™, OrthoWedge®, SlimLine™

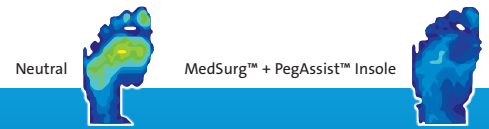

- > for selective off-loading and acute treatments
- > separately available for orthopaedic devices shown at left
- > right and left fit

The PegAssist™ Insole is built around a multi-purpose removable peg insole that effectively off-loads the plantar aspect of the foot. Easy removal of individual pegs allows selective off-loading and the Peg Support Board keeps the offloaded area intact. This insole is constructed of an 18 mm combination of Plastazote®, Poron® and EVA to provide the perfect combination of shock absorption and moldability.

WCS™

Wound Care Shoe System



exchangeable insole

WCS™ Insoles

For further treatment and if lesions are changing, the four layers can be ordered separately.

- > postoperative for diabetic patients
- > for off-loading of ulcerations and wounds
- > 7 Sizes: XS, S, S / M, M, L, XL, XXL

The WCS™ is recommended for the treatment of open and closed ulcerations and other conditions of the foot in which it is desirable to redistribute weight away from specific areas. The circumferential counter of the shoe forms a deep pocket in the sole to permit the use of a variety of insoles under the ulceration or area of pressure. The insole material is placed below the level of the top of the circumferential counter, providing greater stability for the foot by preventing the layered insoles from shifting within the shoe.

WCS™ Closed-Toe

Wound Care Shoe System



exchangeable insole

WCS™ Insoles

For further treatment and if lesions are changing, the four layers can be ordered separately.



- > for treatment of diabetic foot syndrome
- > for healing of open wounds and minor ulcerations
- > for pressure relief
- > wide toe box for any remaining bandages

> 7 Sizes: XS, S, S / M, M, L, XL, XXL

The WCS™ Closed-Toe system can be customized individually for off-loading. The proven quadruple layered and removable insoles enable a selective creation of plantar off-loading zones. The high heel counter provides rear foot stability and the toe upper is made of a soft, elastic, seamless and breathable material for a comfortable fitting. The extra space in the toe area accommodates light bandages and deformities. The WCS™ Closed-Toe incorporates a rocker bottom sole to provide a physiologically correct gait. For indoor and outdoor.

HeelWedge™

Off-loading Shoe



- > for the treatment of plantar fasciitis, ulcerations, infections, trauma and following surgery of the soft tissue or bony structure of the heel

> 5 Sizes: XS, S, M, L, XL
> right and left fit

The HeelWedge™ is clinically proven to off-load pressure by over 26 % from the heel by shifting weight to the mid and forefoot to promote faster healing after surgery, trauma or when wounds or ulcerations are present on the heel. The square-toe design offers more space for bandage material. The modified and removeable forefoot closure provides better left / right fit. The removable EVA-insole can be modified as needed.

compatible with:



TwinShoe
for height adjustment



PegAssist™ Insole
for selective off- loading

exchangeable insole

OrthoWedge®

Off-loading Shoe



- > use after surgeries to correct hallux valgus, hammer toes /
Tailor's bunions
- > treatment of diabetic ulcerations
- > 5 Sizes: XS, S, M, L, XL
- > right and left fit

compatible with:



TwinShoe
for height adjustment



PegAssist™ Insole
for selective off- loading

exchangeable
insole

The OrthoWedge® Off-loading Shoe is specifically designed to protect the forefoot by removing most of the pressure from the metatarsal head and digits during the weight-bearing portion of the gait cycle. For modification the insole can be taken out. The patented construction of the wedge obtains an off-loading mechanism (TPR-sole). The weight is noticeably reduced at the forefoot. Velcro straps make opening easy and allow a comfortable fitting.

TwinShoe

Balance Shoe



- > for a quick mobilization after foot surgery
- > height adjusting balance shoe matches all orthopaedic
devices with a thick sole
- > 5 Sizes: XS, S, M, L, XL
- > Height of the sole: 36 mm
- > right and left fit

compatible with:

**HeelWedge™, OrthoWedge®,
Body Armor® Walker II**

The TwinShoe is a matching and height adjusting shoe supplementing our classic product line (listed bottom left) and used specifically for clinical care support. No more time-intensive applying and removing of height adjustments. DARCO TwinShoe warrants a safe and physiologically correct gait.

MedSurg™


Postoperative Shoe



compatible with:

PegAssist™ Insole
for selective off-loading

exchangeable
insole

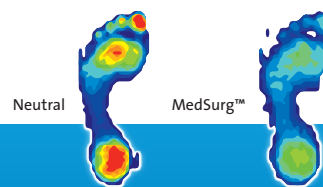

> used postoperatively for stabilization and off-loading

- > 3 Women-Sizes: WS, WM, WL
- > 4 Men-Sizes: MS, MM, ML, MXL
- > right and left fit

The MedSurg™ shoe with its famous semi-rigid sole is built around our patented MetaShank™ forefoot support, that provides more rigid control and better protection after osteotomies. Our breathable, mesh upper conforms to bony abnormalities of the foot and has a reinforced heel counter and padded heel collar to stabilize the rear foot and reduce slippage. Extra-long straps accommodate even the bulkiest dressings. The MedSurg™ shoe is available in 7 different sizes. Helping you get the best possible fitting post-operative shoe.

GentleStep™

Diabetic Shoe



- > ready-to-use protective shoe for diabetics
- > mid rocker sole permits easy ambulation
- > fully breathable Lycra® toe elastic upper reduces weight while increasing comfort
- > firm heel counter ensures stable and comfortable fit
- > machine washable

> 19 Sizes: available in nineteen sizes and two widths

The GentleStep™ Shoe is an extra-depth shoe designed specifically for the needs of certain patients with diabetes. The upper is constructed of a light-weight Lycra® to provide extra comfort for patients with forefoot deformities or lesions.

exchangeable
insole



GentleStep™ Insole

Off-loading insole



compatible with:

GentleStep™

> for selective off-loading and acute treatments

> separately available for orthopaedic devices shown at left

> fits to existing GentleStep™ sizes

> sold in pairs

The GentleStep™ Insole is built around a multi-purpose removable peg insole that effectively off-loads the plantar aspect of the foot. Easy removal of individual pegs allows selective off-loading and the Peg Support Board keeps the off-loaded area intact. This insole is constructed of an 10 mm combination of Poron and EVA to provide the perfect combination of shock absorption and moldability.

WCS™/-Light Insoles

4-Layer Insole System



compatible with:

WCS™ Wound Care Shoe and
WCS™ Closed Toe Wound Care
Shoe

- > treatment and management of diabetic foot syndrome
- > treatment for healing of open wounds, ulcers and ulcerations
- > postoperative phase after toe, forefoot and midfoot surgical procedures
- > select customizable pressure relief

The WCS™ insole provides plantar pressure relief through the easy task of adapting the Four-Layer-Insole-System. The four-layer system consists of EVA, PE and Poron®. Each differ through varying grades stiffness and elastic recovery.



Plantar Fasciitis Products

Pneumatic Night Splint

Posterior Night Splint



- > prophylaxis of drop foot
- > tendinitis in heel area
- > for immobilization
- > 4 Women-Sizes: S, M, L, XL
- > 3 Men-Sizes: M, L, XL
- > right and left fit

The Pneumatic Night Splint features a fully-integrated inflation / deflation system that can be easily operated with one hand. The internal air bladder is infinitely adjustable so every wearer can find the degree of dorsiflexion that works best for them. The natural cotton terry liner absorbs moisture and ensures a cool, comfortable fit and a non-skid sole for minimal ambulation when deflated.

Body Armor® Night Splint

Dorsal Night Splint




inclusive:

Toe Plate

The toe plate effectively places the toes in a dorsiflexion, promoting the stretching of the plantar fascia.

- > shortening of the plantar fasciitis
- > heel spurs and metatarsalgia
- > achilles tendonitis and tendonosis
- > talipes equinus and ankle contracture
- > 1 Size: Universal
- > right and left fit

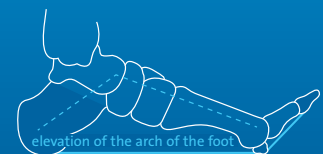

The dorsal night splint stabilizes the foot and lower leg position at an angle of 90° in order to avoid equinus deformity of the foot and to consequently prevent achilles tendon shortening. The toe loop elevates the hallux (dorsal flexion) resulting in a sustained stretch to the plantar fascia. The combination of stretching to the fascia and keeping the position by means of the splint enables various indications to counteract. The adjustable straps are quick and easy to use providing perfect fitting.

DCS™ Plantar Fasciitis Sleeve

Plantar Fasciitis Sleeve



- > wave compression provides a consistent and sustained stretch to the plantar fascia
- > for Plantar Fasciitis, Heel Spurs, Arch pain, Arch weakness, Edema
- > 3 Colors: blue, black, tan
- > 4 Sizes: S, M, L, XL
- > right and left fit

The DARCO DCS™ Plantar Fasciitis Sleeve is woven to be anatomically correct. Wave compression technology on the plantar surface of the foot provides a sustained stretch on the plantar fascia. The DCS™ PF gives relief from the discomfort associated with plantar fasciitis by promoting circulation, reducing edema and foot fatigue, and relieving morning foot pain.

ProMotion® Plus

Orthotic Insoles



- > for Over-pronation, Plantar fasciitis, Heel pain, Foot and arch pain, Ankle pain, Knee pain
- > 9 Women-Sizes, 9 Men-Sizes

The DARCO ProMotion® Plus orthotic insole features a full length PORON® Cushioning layer which provides greater comfort and doesn't bottom out long term.



Ankle Supports

Body Armor® Vario

Ankle Support



inclusive:

Heel posting disks

3 Angles of inclination 2°, 4°, 6°

Arch height adjustments

3 Sizes S, M, L



- > for peroneal injuries
- > for lateral ankle ligament sprains
- > for deltoid sprains
- > for posterior tibial tendon injuries
- > after flatfoot reconstruction
- > after lapidus fusion

- > 1 Size: Universal
- > for right and left side available

A DARCO advisory board of certified foot and ankle specialist developed the brace as a treatment option for patients with ankle instabilities along with alignment issues (Varus /Valgus), i.e. flat foot caused by PTTD., as well as for supportive ankle /arch positions. The Body Armor® Vario brace is a general use ankle brace that adds an element of customization not seen in other »off-the-shelf« products.

Web™ Ankle Brace

Ankle Support



Removable medial and lateral rigid stays.

- > for prevention during sports or in everyday life as well as a supportive measure in minor acute ankle injuries

- > 5 Sizes: XS, S, M, L, XL
- > right and left fit

The Web™ Ankle Brace is a true convertible ankle brace. The bungee closure provides superior compression and allows for adjustments on the run. Removable medial and lateral stays are rigid enough to provide controlled restriction of inversion /eversion. Durable ballistic nylon shell is breathable and conforms to the foot for a comfortable fit. The brace is constructed of a durable vinyl shell – easy to get in and thin enough to wear with athletic shoes or casual dress shoes.



Pediatric Products

Pediatric MedSurg™

Postoperative Shoe



- > used postoperatively for stabilization and off-loading
- > 1 Size
- > right and left fit

The Pediatric MedSurg™ is the smaller version of our MedSurg™ postoperative shoe for adults. The child's healing shoe contains all the beneficial features found in the adult MedSurg™ but sized for smaller feet. The quality and construction characteristics remain the same.

Pediatric Original MedSurg™

Postoperative Shoe



- > used postoperatively for stabilization and off-loading
- > 1 Size
- > right and left fit

The DARCO medical-surgical shoe is the shoe that started it all. Our Original MedSurg™ shoe was invented over 30 years ago and has since become known as the industry standard. Often referred to generically as »the DARCO shoe«, it is anything but generic and its long list of features and excellent quality make it one of the best values in the business.

Pediatric SlimLine™

Cast Boot



- > used over a fiberglass- or plaster cast
- > 3 Sizes: S, M, L
- > right and left fit

The Pediatric SlimLine™ Cast Boot offers protection of the cast against dirt and exposure. All features and characteristics of the standard model for adults are included in this product as well.







Order Information

Sizes and Measurements


MedSurg™, Original MedSurg™, Softie™ Shoe (Women)

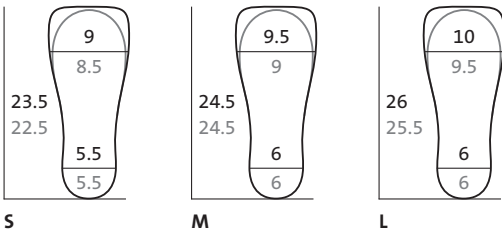


MedSurg™, Original MedSurg™, Softie™ Shoe (Men)

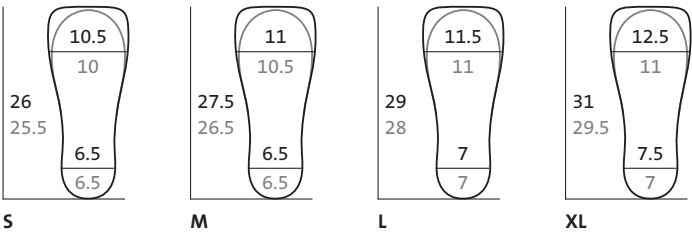


APB™, HeelWedge™, OrthoWedge®, TwinShoe, SlimLine™

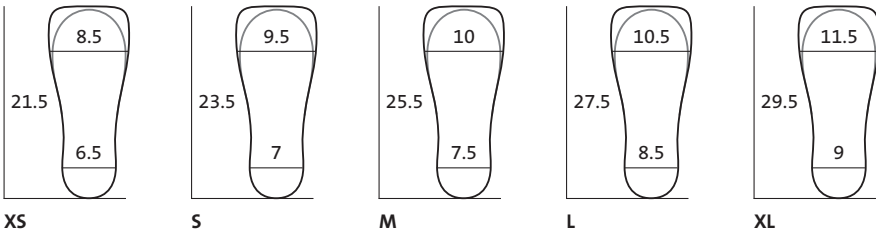


Body Armor® Pro Term

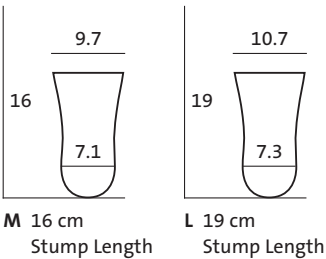


FXPro™ Walker, AirTraveler™

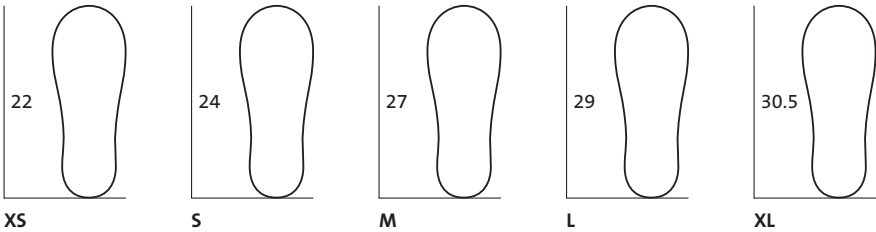

Body Armor® Walker II

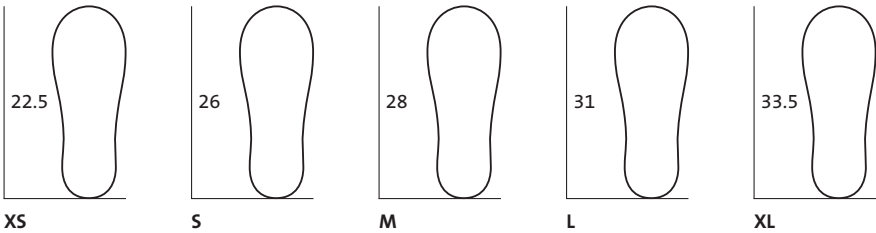


WCS™, WCS™ Closed-Toe

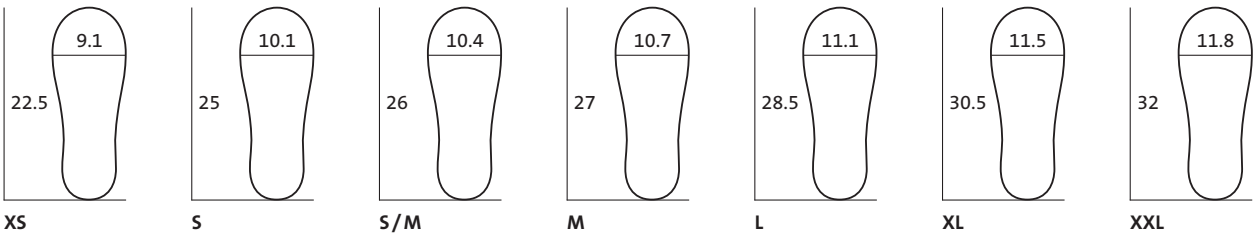

Sizes and Measurements Kids




Pediatric MedSurg™, Pediatric Original MedSurg™

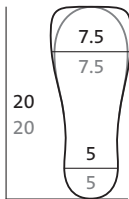


1 Size


Pediatric SlimLine™

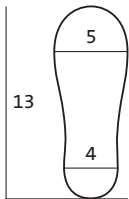

S

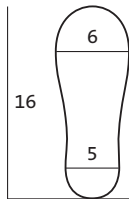

M

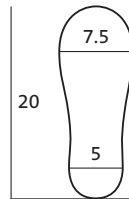

L

All dimensions are shown in cm and relate to the insoles.

Attention

Medical and/or orthopaedic and/or therapeutic footwear needs to be applied by an experienced HCP (Health Care Professional). The shoe sizes mentioned in this catalog are for orientation only. Wound dressings, deformities, presense of implants, healing process, amputations and many other indications might require bigger or smaller shoe sizes as the patient is used to wear with his regular street footwear.

»Experience is required to finally establish that the footwear is the most appropriate fitting and size for the foot. Whatever size scale is printed on the inside of the footwear, or whatever size is measured on a gauge, there is no substitute for the valuable experience of an expert in shoe fitting who can ensure that the shoe and the foot complement each other in size, form and function.«

Tyrrell, W.; Carter, G. (2009) »Therapeutic Footwear – A comprehensive Guide«, Churchill Livingstone, p. 84

Article List

Item-No. Colour Size (UK) Women Men

Postoperative Care



MedSurg™

Women:				
MQW1B	Black	S	2.5–3.5	
MQW2B	Black	M	4.0–5.0	
MQW3B	Black	L	5.5–7.5	
Men:				
MQM1B	Black	S		6.0–8.0
MQM2B	Black	M		8.5–10.0
MQM3B	Black	L		10.5–12.5
MQM4B	Black	XL		13.0–14.0



Original MedSurg™

Women:				
MSW1N	Blue	S	1.0–2.5	
MSW2N	Blue	M	3.0–5.0	
MSW3N	Blue	L	5.5–6.5	
Men:				
MSM1N	Blue	S		5.5–7.0
MSM2N	Blue	M		7.5–8.5
MSM3N	Blue	L		9.0–11.0
MSM4N	Blue	XL		11.5–13.5



Softie™ Shoe

Women:				
STW1B	Black	S	1.0–2.5	
STW2B	Black	M	3.0–5.0	
STW3B	Black	L	5.5–6.5	
Men:				
STM1B	Black	S		5.5–7.0
STM2B	Black	M		7.5–8.5
STM3B	Black	L		9.0–11.0
STM4B	Black	XL		11.5–13.5



APB™

APQ0B	Black	XS	1.0–1.5	1.0–2.0
APQ1B	Black	S	2.0–3.5	2.5–4.0
APQ2B	Black	M	4.0–6.5	4.5–7.0
APQ3B	Black	L	7.0–9.5	7.5–10.0
APQ4B	Black	XL	10.0–13.0	10.5–13.0



HeelWedge™

HQ0B	Black	XS	1.0–1.5	1.0–2.0
HQ1B	Black	S	2.0–3.5	2.5–4.0
HQ2B	Black	M	4.0–6.5	4.5–7.0
HQ3B	Black	L	7.0–9.5	7.5–10.0
HQ4B	Black	XL	10.0–13.0	10.5–13.0



OrthoWedge®

OQ0B	Black	XS	1.0–1.5	1.0–2.0
OQ1B	Black	S	2.0–3.5	2.5–4.0
OQ2B	Black	M	4.0–6.5	4.5–7.0
OQ3B	Black	L	7.0–9.5	7.5–10.0
OQ4B	Black	XL	10.0–13.0	10.5–13.0



TwinShoe

TS0B	Black	XS	1.0–1.5	1.0–2.0
TS1B	Black	S	2.0–3.5	2.5–4.0
TS2B	Black	M	4.0–6.5	4.5–7.0
TS3B	Black	L	7.0–9.5	7.5–10.0
TS4B	Black	XL	10.0–13.0	10.5–13.0

Height of the sole: 36 mm



TAS™ (Toe Alignment Splint)

TAS	White	Universal
-----	-------	-----------



GTS™

Women:				
Right	Black			
GTS1R		S	5 – 7.5	
GTS2R		L	8 – 11	
Left	Black			
GTS1L		S	5 – 7.5	
GTS2L		L	8 – 11	
Men:				
Right	Black			
GTS1R		S		3.5 – 6
GTS2R		L		6.5 – 9.5
Left	Black			
GTS1L		S		3.5 – 6
GTS2L		L		6.5 – 9.5

Item-No.	Colour	Size (UK)	Women	Men
----------	--------	-----------	-------	-----

Trauma Care



Body Armor® Cast Shoe

BACS0	Dark Grey	XS	1.5–2.5	2.0–3.0
BACS1	Dark Grey	S	3.0–6.5	3.5–7.0
BACS2	Dark Grey	M	7.0–10.5	7.5–11.0
BACS3	Dark Grey	L	11.0–13.5	11.5–13.5
BACS4	Dark Grey	XL	> 13.5	> 13.5



SlimLine™

SLQ0B	Black	XS	1.0–1.5	1.0–2.0
SLQ1B	Black	S	2.0–3.5	2.5–4.0
SLQ2B	Black	M	4.0–6.5	4.5–7.0
SLQ3B	Black	L	7.0–9.5	7.5–10.0
SLQ4B	Black	XL	10.0–13.0	10.5–13.0



FXPro™ Walker

high version:

FX0	Black	XS	1.5–2.0	2.0–2.5
FX1	Black	S	2.5–4.0	3.0–4.5
FX2	Black	M	4.5–8.0	5.0–9.0
FX3	Black	L	8.0–11.5	9.5–12.5
FX4	Black	XL	11.5–13.5	13.0–14.5

low version:

FXS0	Black	XS	1.5–2.0	2.0–2.5
FXS1	Black	S	2.5–4.0	3.0–4.5
FXS2	Black	M	4.5–8.0	5.0–9.0
FXS3	Black	L	8.0–11.5	9.5–12.5
FXS4	Black	XL	11.5–13.5	13.0–14.5



AirTraveler™

high version:

AT0	Black	XS	1.5–2.0	2.0–2.5
AT1	Black	S	2.5–4.0	3.0–4.5
AT2	Black	M	4.5–8.0	5.0–9.0
AT3	Black	L	8.0–11.5	9.5–12.5
AT4	Black	XL	11.5–13.5	13.0–14.5

low version:

ATS0	Black	XS	1.5–2.0	2.0–2.5
ATS1	Black	S	2.5–4.0	3.0–4.5
ATS2	Black	M	4.5–8.0	5.0–9.0
ATS3	Black	L	8.0–11.5	9.5–12.5
ATS4	Black	XL	11.5–13.5	13.0–14.5

Item-No.	Colour	Size (UK)	Women	Men
----------	--------	-----------	-------	-----



Body Armor® Walker II

high version:

BAT0	Black	XS	1.5–2.0	1.5–2.5
BAT1	Black	S	2.5–7.5	3.0–8.0
BAT2	Black	M	8.0–10.5	8.5–11.0
BAT3	Black	L	11.0–13.5	11.5–14.0
BAT4	Black	XL	> 13.5	> 14.0

low version:

BATS0	Black	XS	1.5–2.0	1.5–2.5
BATS1	Black	S	2.5–7.5	3.0–8.0
BATS2	Black	M	8.0–10.5	8.5–11.0
BATS3	Black	L	11.0–13.5	11.5–14.0
BATS4	Black	XL	> 13.5	> 14.0



Body Armor® ProTerm

BAPT2	Black	M	16.0 cm	16.0 cm
BAPT3	Black	L	19.0 cm	19.0 cm

Off-loading and Wound Care



PegAssist™ Insole (PTQ Series for MedSurg™)

Women:

PTQW1	S	2.5–3.5
PTQW2	M	4.0–5.0
PTQW3	L	5.5–7.5

Men:

PTQM1	S	6.0–8.0
PTQM2	M	8.5–10.0
PTQM3	L	10.5–12.5
PTQM4	XL	13.0–14.0



PegAssist™ Insole (PQ Series for APB™, HeelWedge™, OrthoWedge®, SlimLine™)


PQ0	XS	1.0–1.5	1.0–2.0
PQ1	S	2.0–3.5	2.5–4.0
PQ2	M	4.0–6.5	4.5–7.0
PQ3	L	7.0–9.5	7.5–10.0
PQ4	XL	10.0–13.0	10.5–13.0

Item-No. Colour Size (UK) Women Men



WCS™

WCS0B	Black	XS	1.5–2.0	2.0–2.5
WCS1B	Black	S	2.5–5.5	3.0–6.0
WCS7B	Black	S/M	6.0–7.5	6.5–8.0
WCS2B	Black	M	8.0–8.5	8.5–9.5
WCS3B	Black	L	9.0–10.5	10.0–11.5
WCS4B	Black	XL	11.0–13.5	12.0–13.5
WCS5B	Black	XXL	> 13.5	> 13.5

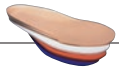


+ Insoles (for WCS™, 4 Pieces)

WCI0	XS	1.5–2.0	2.0–2.5
WCI1	S	2.5–5.5	3.0–6.0
WCI7	S/M	6.0–7.5	6.5–8.0
WCI2	M	8.0–8.5	8.5–9.5
WCI3	L	9.0–10.5	10.0–11.5
WCI4	XL	11.0–13.5	12.0–13.5
WCI5	XXL	> 13.5	> 13.5



WCS™ Closed-Toe

WCC0B	Black	XS	1.5–2.0	2.0–2.5
WCC1B	Black	S	2.5–5.5	3.0–6.0
WCC7B	Black	S/M	6.0–7.5	6.5–8.0
WCC2B	Black	M	8.0–8.5	8.5–9.5
WCC3B	Black	L	9.0–10.5	10.0–11.5
WCC4B	Black	XL	11.0–13.5	12.0–13.5
WCC5B	Black	XXL	> 13.5	> 13.5

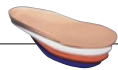

+ Insoles (for WCS™ Closed-Toe, 4 Pieces)

WCI0	XS	1.5–2.0	2.0–2.5
WCI1	S	2.5–5.5	3.0–6.0
WCI7	S/M	6.0–7.5	6.5–8.0
WCI2	M	8.0–8.5	8.5–9.5
WCI3	L	9.0–10.5	10.0–11.5
WCI4	XL	11.0–13.5	12.0–13.5
WCI5	XXL	> 13.5	> 13.5



HeelWedge™

HQ0B	Black	XS	1.0–1.5	1.0–2.0
HQ1B	Black	S	2.0–3.5	2.5–4.0
HQ2B	Black	M	4.0–6.5	4.5–7.0
HQ3B	Black	L	7.0–9.5	7.5–10.0
HQ4B	Black	XL	10.0–13.0	10.5–13.0

Item-No. Colour Size (UK) Women Men



OrthoWedge®

OQ0B	Black	XS	1.0–1.5	1.0–2.0
OQ1B	Black	S	2.0–3.5	2.5–4.0
OQ2B	Black	M	4.0–6.5	4.5–7.0
OQ3B	Black	L	7.0–9.5	7.5–10.0
OQ4B	Black	XL	10.0–13.0	10.5–13.0



TwinShoe

TS0B	Black	XS	1.0–1.5	1.0–2.0
TS1B	Black	S	2.0–3.5	2.5–4.0
TS2B	Black	M	4.0–6.5	4.5–7.0
TS3B	Black	L	7.0–9.5	7.5–10.0
TS4B	Black	XL	10.0–13.0	10.5–13.0



MedSurg™

Women:

MQW1B	Black	S	2.5–3.5
MQW2B	Black	M	4.0–5.0
MQW3B	Black	L	5.5–7.5

Men:

MQM1B	Black	S	6.0–8.0
MQM2B	Black	M	8.5–10.0
MQM3B	Black	L	10.5–12.5
MQM4B	Black	XL	13.0–14.0



GentleStep™

Wide (D & E):

	Women	Men
GS070WB Black	4.5	5.0
GS075WB Black	5.0	5.5
GS080WB Black	5.5	6.0
GS085WB Black	6.0	6.5
GS090WB Black	6.5	7.0
GS095WB Black	7.0	7.5
GS100WB Black	7.5	8.0
GS105WB Black	8.0	8.5
GS110WB Black	8.5	9.0
GS115WB Black	9.0	9.5
GS120WB Black	9.5	10.0
GS125WB Black	10.0	10.5
GS130WB Black	10.5	11.0
GS135WB Black	11.0	11.5
GS140WB Black	11.5	12.0
GS145WB Black	12.0	12.5
GS150WB Black	12.5	13.0
GS155WB Black	13.0	13.5
GS160WB Black	13.5	14.0

Extra-Wide (EE & EEE):

	Women	Men
GS070XB Black	4.5	5.0
GS075XB Black	5.0	5.5
GS080XB Black	5.5	6.0
GS085XB Black	6.0	6.5

Item-No.	Colour	Size (UK) Women	Men
----------	--------	-----------------	-----

GS090XB	Black	6.5	7.0
GS095XB	Black	7.0	7.5
GS100XB	Black	7.5	8.0
GS105XB	Black	8.0	8.5
GS110XB	Black	8.5	9.0
GS115XB	Black	9.0	9.5
GS120XB	Black	9.5	10.0
GS125XB	Black	10.0	10.5
GS130XB	Black	10.5	11.0
GS135XB	Black	11.0	11.5
GS140XB	Black	11.5	12.0
GS145XB	Black	12.0	12.5
GS150XB	Black	12.5	13.0
GS155XB	Black	13.0	13.5
GS160XB	Black	13.5	14.0

Item-No.	Colour	Size (UK) Women	Men
----------	--------	-----------------	-----

ProMotion® Plus

PMP03	Blue	3	5 – 5.5	–
PMP04	Blue	4	6 – 6.5	–
PMP05	Blue	5	7 – 7.5	5 – 5.5
PMP06	Blue	6	8 – 8.5	6 – 6.5
PMP07	Blue	7	9 – 9.5	7 – 7.5
PMP08	Blue	8	10 – 10.5	8 – 8.5
PMP09	Blue	9	11 – 11.5	9 – 9.5
PMP10	Blue	10	12 – 12.5	10 – 10.5
PMP11	Blue	11	13 – 13.5	11 – 11.5
PMP12	Blue	12	–	12 – 13.5
PMP13	Blue	14	–	14 – 14.5

GentleStep™ Insole

	for GentleStep™ Size	
	Women	Men
GSI 36	up to 7.0	up to 5.5
GSI 37	7.5	6.0
GSI 38	8.0 – 8.5	6.5 – 7.0
GSI 39	9.0 – 9.5	7.5 – 8.0
GSI 40	10.0	8.5
GSI 41	10.5	9.0
GSI 42	11.0 – 11.5	9.5 – 10.0
GSI 43	12.0	10.5
GSI 44	12.5	11.0
GSI 45	13.0	11.5
GSI 46	13.5	12.0

Ankle Supports

Body Armor® Vario

left:		
BAV2-L	Black	Universal
right:		
BAV2-R	Black	Universal

Web™ Ankle Brace

WB0	Black	XS	2.5–3.5	3.0–4.0
WB1	Black	S	4.0–6.5	4.5–7.0
WB2	Black	M	7.0–8.0	7.5–8.5
WB3	Black	L	9.0–10.5	9.5–11.0
WB4	Black	XL	11.0–13.0	11.5–13.5

Plantar Fasciitis Products

Pneumatic Night Splint

NS1B	Black	S	1.5–2.5	2.0–3.0
NS2B	Black	M	3.0–5.5	3.5–6.0
NS3B	Black	L	6.0–8.0	6.5–9.5
NS4B	Black	XL	8.5–14.0	10.0–14.5

Body Armor® Night Splint

BADS	White-Blue	Universal
------	------------	-----------

DCS™ Plantar Fasciitis Sleeve

	Arch Circumference (inch)	
DCS-PF1	White-Blue	S 5.0 – 8.0
DCS-PF2	White-Blue	M 7.0 – 10.0
DCS-PF3	White-Blue	L 9.0 – 13.0
DCS-PF4	White-Blue	XL 11.0 – 15.0

Pediatric Products

Pediatric MedSurg™

MQPB	Black	1 Size	< 12.5
------	-------	--------	--------

Pediatric Original MedSurg™

MSPN	Blue	1 Size	< 12.0
------	------	--------	--------

Pediatric SlimLine™

SLO1NP	Blue	S	< 6.0
SLO2NP	Blue	M	6.5–9.0
SLO3NP	Blue	L	9.5–12.5

Legal Notice

DARCO does not practise in the medical field, nor does DARCO make recommendations for specific application to patients. The physician in attendance is responsible for deciding on the appropriate medicine and for the carrying out of the appropriate treatment in every individual patient case. DARCO assumes no responsibility for this.

Publisher and responsible for contents

DARCO Medical India Pvt. Ltd.
32A1B1/A1, First Floor, Vaibhav Complex
Opp to District Court, P. B. Road
Dharwad 580 008, Karnataka, India

Telephone +91. 836. 2 74 74 46
Mobile +91. 94 83 81 44 22
info@darco.in, www.darco.in

Last update

March 2019

© DARCO 2019

Concept and design

westcoastmedia
www.westcoastmedia.de

Product photography

DARCO (Europe) GmbH